



Hyber Hybridoma Cell Serum-Free Media

A serum-free medium especially developed for high-density suspension culture of hybridoma cells

Hyber hybridoma cell medium is a serum-free medium independently developed by Shanghai BioEngine Sci-Tech Co., Ltd., and it is an ideal choice for culturing hybridoma cells from laboratory to large-scale suspension. It enables rapid growth of hybridoma cells and supports efficient expression of antibodies in hybridoma cells while providing a guarantee for the safety and reliability of the production process.

Features

- A completely serum-free culture system
- No ingredients of animal-derived, genetically modified plant origin or raw materials with BSE origin
- No antibiotics, organic solvents, or preservatives
- Containing hydrolysate, L-glutamine, hypoxanthine, and thymine nucleoside
- Supporting high-density suspension culture and protein expression of hybridoma cells



Hyber Hybridoma Cell Serum-free Media

Advantages

- Animal-derived component-free; TSE/BSE statement available on demand;
- Distinctive culture results proven in numerous studies;
- Optional powder media for use in large-scale manufacturing with easy preparation procedures;
- Powder media capable of a single batch size of 100,000 L;
- Excellent inter-batch consistency (CPK*>1.33);
- Full traceability by EU-certified ISO13485:2016 Quality Management System.

*CPK: Process Capability Index; a CPK > 1.33 indicates good process control and small inter-batch difference of products.


Ordering Information

Product Name	Cat. No.	Form	Size	Package	Other	
Hyber-B100S Hybridoma Cell Serum-free Medium	EXP0111201	Powder	100L	Bag	[+]L-GIn [-]NaHCO₃	
	EXP0111202	Powder	10L	Bag	[+]L-GIn [-]NaHCO₃	
	EXP0111203	Powder	5L	Bag	[+]L-GIn [-]NaHCO ₃	
Hyber-B100 Hybridoma Cell Serum-free Medium	EXP0101402	Liquid	1L	Bottle	[+]L-GIn [+]NaHCO₃	
Hyber-F100S Hybridoma Cell Serum-free Feed Medium	EXP0111301	Powder	20L	Bag	[+]L-GIn [-]NaHCO₃	
	EXP0111302	Powder	2L	Bag	[+]L-GIn [-]NaHCO₃	
Hyber-F100 Hybridoma Cell Serum-free Feed Medium	EXP0100802	Liquid	250ml	Bottle	[+]L-GIn [+]NaHCO₃	

Performance


Suspension acclimation of adherent hybridoma cells

The use of Hyber-B100S with gradual reduction of the serum content during continued acclimation of hybridoma cells allowed acclimation of serum-adherent cells into serum-free suspension cells. Suspension cells were individually dispersed with a uniform cell size and showed a stable specific growth rate over the passages.

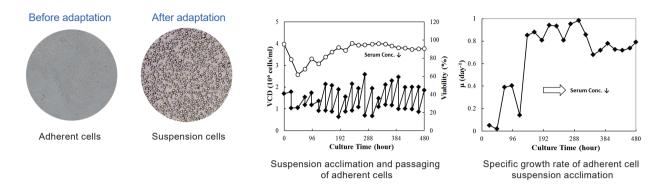

Satch culture of suspension hybridoma cells

Hyber-B100S is used for batch culture of hybridoma cells with a maximum viable cell density of 8×10⁶ cells/ml and is more effective when used with Hyber-F100S.

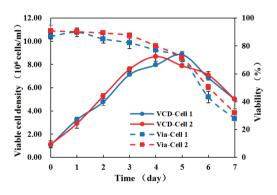

Fed-batch usage strategy

Instructions for use of Hyber-F100S:

Culture time	C0	C1	C2	C3	C4	C5	C6	C7
Ratio of <i>Hyber-F100S</i> to C0 initial medium v/v (ml/ml)	0	3%	3%	3%	3%	3%	3%	0

Tel ·

30 years of ingenuity on creating a novel drive for cell culture



Add.:

Floor 5, Building 1, Lane 720, Cailun Road, Zhangjiang Hi-Tech Park, Pudong, Shanghai, PRC

www.bio-engine.com.cn (86)21-68582660

Web .:

E-mail: ${\bf extsf{matrix}}$ marketing@bio-engine.com.cn